


EYE & TEMPLE MASSAGER

Instruction Manual


Model: EYE03

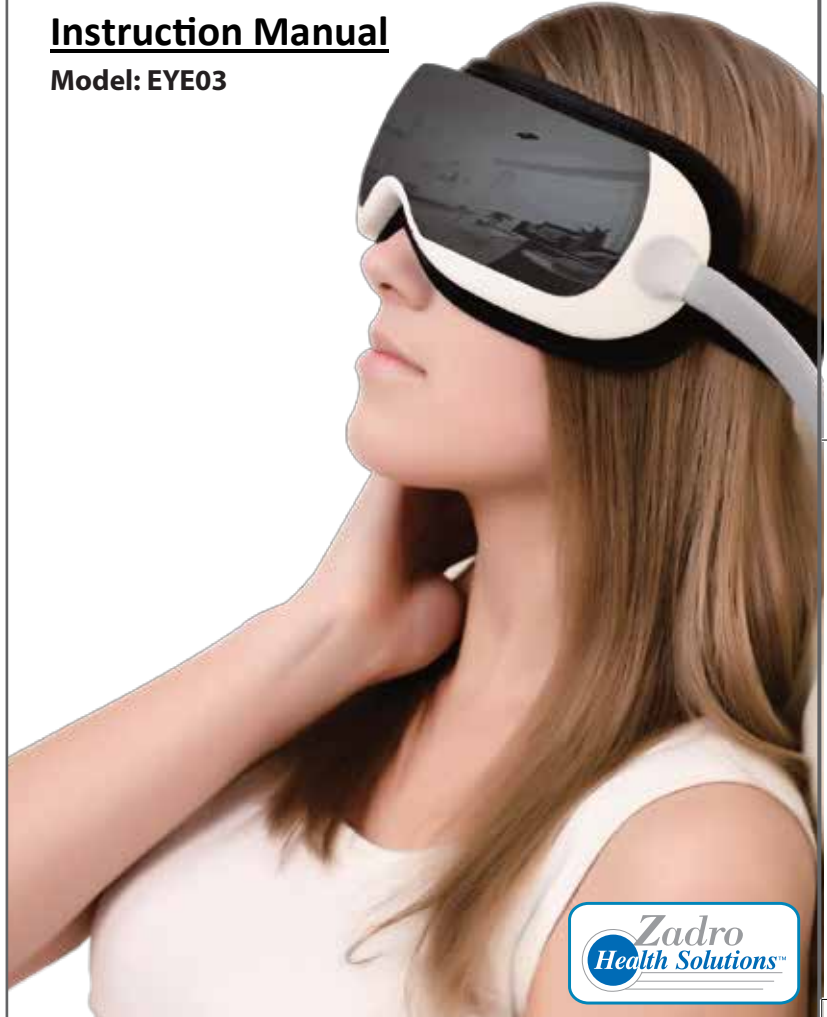

Zadro Health Solutions, Inc.
Innovative Ideas That Make Life Easier™
Huntington Beach, CA 92649
Phone: (714) 892-9200 • 800-468-4348
Fax: (714) 892-7477
E-mail: ZadroHS@ZadroHS.com
Website: www.ZadroHS.com

Rev.04232014



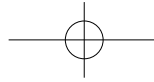

Thank you for choosing Zadro Health Solutions Eye & Temple Massager. Please visit us online to leave comments, reviews, & see other great products from www.Zadrolnc.com

Please read the Instruction Manual carefully before use. Retain this Instruction Manual for future reference.

Table of Contents

Introduction and Benefits	2
Product Description	2
Product Diagram	3
Product Functions	4
Safety Information	5
Product Setup	6
Operation Method	7
Care and Handling	8
Troubleshooting	8
Specifications	9
Warranty	10

Introduction and Benefits

Using the Zadro Health Solutions Eye and Temple Massager will help you alleviate eye and facial strain from everyday stress allowing you to relax and be more at ease.

The Eye and Temple Massager uses eleven different vibration settings to help relieve tired eyes from the stress of computer usage, reading, driving, and more. Massaging the acupressure points around the eyes has been shown to relieve eye and facial strain. This product also has two heat settings. Warming the area around the eyes will help relax the ciliary muscles. Having a warm, soft massage will relax your eyes helping them feel more refreshed and relaxed. In addition, there are five different types of air pressure generated by a hi-tech IC chip which helps improve your microcirculation in blood vessels around the eyes. Using the Eye and Temple Massager is a smart way to relax.


Eye and temple massage relaxes your muscles helping to relieve insomnia, headaches, and other tension and stress related issues.

An additional bonus with this product are the included earphones allowing you to listen to the provided music or upload your own MP3 files to the unit.


Product Description

- For a Customizable, Gentle Eye Massage there are 11 Types Of Vibrations, Air Pressure, and Heat Settings
- The Eye and Temple Massager Relaxes Your Eyes and Relieves Eye Fatigue and Facial Strain.
- Adjustable Head Strap Fits The Eye and Temple Massager Comfortably on the Head.
- Convenient and Portable. Uses Four AA Batteries or AC Adapter.

Product Diagram

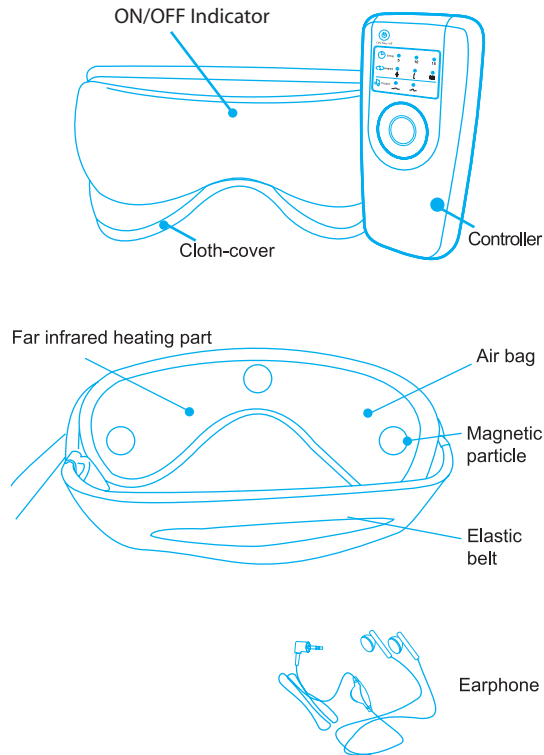

Product Functions

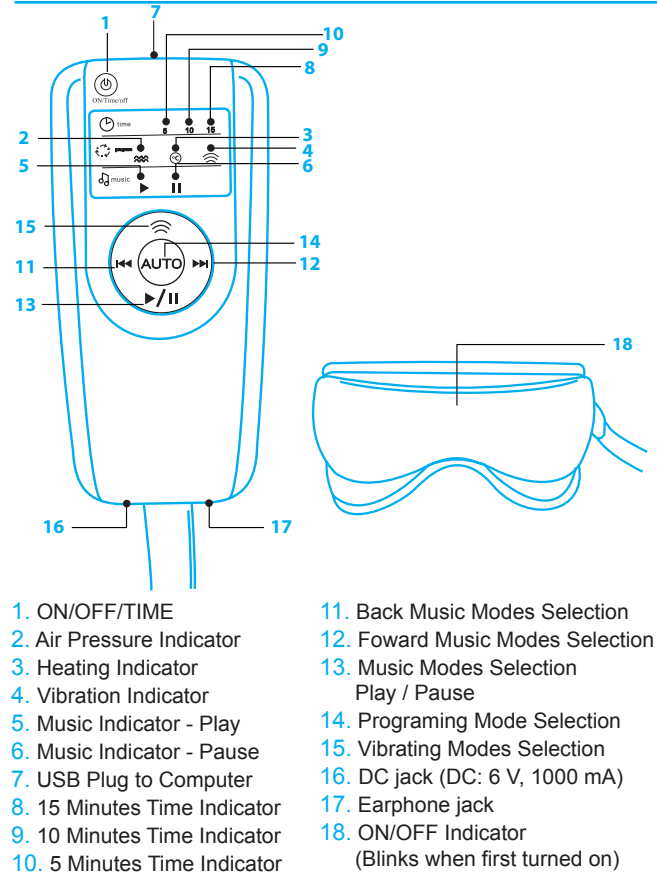

Safety Information

- Carefully read through the operating manual before using the Eye & Temple Massager.
- Do not use this product while bathing or around water. Do not immerse in water or expose to high humidity.
- This product should not be dropped. Please avoid exposing this product to strong sunlight or corrosive gas.
- Eye & Temple Massager covers the eyes so it cannot be used while driving, walking, or in any situation where vision cannot be obstructed.
- Contact lenses must be removed from the eyes before using the Eye & Temple Massager.
- Place the product properly over face and eyes before turning on the product.
- After using the Eye & Temple Massager, keep your eyes closed for 1 to 3 minutes before opening. Then open slowly after use.
- Use alkaline batteries or AC adapter. Do not take out batteries when you are using adapter.
- In case of mechanical malfunction, please contact the manufacturer. The warranty is voided if the product has been dismantled or tampered with.


Do not use this product in the following situations:

- When wearing contact lenses.
- If any medical concerns or when medical conditions contradict use of this product such as glaucoma, retinoschisis, or other eye diseases, eye problems or medical issues, please contact your physician prior to using this product.
- Check eye mask and hose regularly. Do not use this product in the case of a broken or fractured eye mask or cable.

Product Setup

Battery installation


1. Open the battery cover of the controller and insert 4 AA batteries properly.
2. Close the battery cover of the controller.

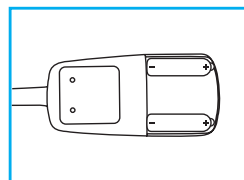

- ⚠ Caution:
Check the direction of the batteries. Install the batteries correctly; check the battery cover and keep it secure.

Adapter installation


1. Connect the adapter with AC power jack properly.
2. Plug the DC PIN into the DC jack of the controller.

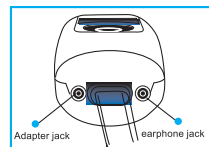


- 1 Connect the earphone with the jack.



- 2 Wear the eye mask.

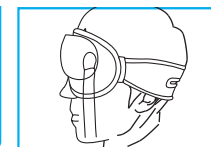

- 3 Adjust the elastic to proper position and fix it around the head.

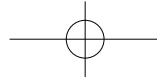

Operation Method

Press On/Off button, the device starts to work



When the device starts, working indicator on the front of mask blinks. Starting indicators showing air pressure, heating, relaxing music, and time for 5 minutes will be on.

Using the Programming mode selection button you can choose different features such as air pressure + heating, vibrating + heating, air pressure + vibrating. Default setting is air pressure massage on temple and eyes, heating is on, timer is set for 5 minutes and the music is on. To customize the settings, follow the charts.


Music:

Music Controls:
Press music button Forward to Skip to next song.
Press music button Back to replay or go back to last song
Press Play/Pause to Stop or Play Music

Music Upload:
Connect USB Mini-B cable to top of remote and computer.
Follow instructions provided on the computer to upload music to EYE03.


Profiles:

Press Programming mode selection button, profile for pressure, heat, and vibration will change

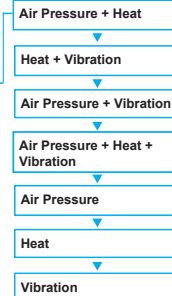

Timer:

Press Time button, Time settings will change

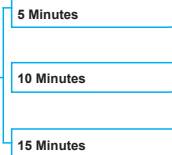

Vibration:

Press Vibration button, Vibration settings will change

High, Medium, Low
Kneading, Circular

Pulsing:
Heavy, Medium, Light

Alternating:
Fast, Slow

Care and Handling

- Do not use this product while bathing or around water. Do not immerse in water or expose to high humidity. Do not wash the Eye & Temple Massager with water. If necessary, clean gently using alcohol.
- This product should not be dropped. Please avoid exposing this product to strong sunlight or corrosive gas.
- Check eye mask and hose regularly. Do not use this product in the case of a broken or fractured eye mask or cable.
- Use alkaline batteries or AC adapter. Do not take out batteries when you are using adapter.

Battery Replacement

1. Use only fresh batteries of the required type and size.
2. Do not mix old and new batteries, different types of batteries, or rechargeable batteries of different capacities.
3. The batteries may explode if mistreated. Do not attempt to recharge or disassemble the batteries.
4. Depleted batteries are to be removed from the item.
5. Dispose of the old batteries properly, and according to your local regulations.
6. Ensure that batteries are kept out of reach of young children.

Troubleshooting

In case of mechanical malfunction, please contact Customer Service at Zadro Health Solutions, Inc. 1 (800) 468-4348 or visit us on the web at www.ZadroHS.com

The warranty is voided if the product has been dismantled or tampered with.



Specifications

Specifications:

Operating Environment:

- Temperature: 41 - 104° F (5 to 40° C)
- Relative Humidity: less than 80%
- Air pressure: 12.5 psi - 15.4 psi (860 hPa to 1060 hPa)
- Intensity of Magnetic Field: greater than 2 kG (200 mT)
- Vibration Frequency: 40 to 100 Hz
- Heating Temperature: 86 - 113° F (30 to 45° C)

Power Supply:

Operating Environment:

AC: 110 V – 240 V, 50 Hz – 60 Hz (adapter optional)

DC: 6 V

Power: <6 W

Weight: >0.5 kg

Controller dimension: 5.5" x 2.8" x 1.6" (14cm x 7cm x 4cm) (L x W x H)

Main Mask dimension: 12" x 4.3" x 4" (30cm x 11cm x 10cm) (L x W x H)

Warranty

90 day limited warranty from date of purchase against defects in materials and workmanship. This warranty does not include damage to the product as a result of misuse or accident.

Should your product not function as expected, please contact our Customer Service department at customerservice@zadroinc.com. When return authorization is obtained, please be sure to package the item securely and return it, prepaid/insured, along with proof of purchase to: Zadro Health Solutions, Inc. Attn: Customer Service, 5422 Argosy Ave. Huntington Beach, CA 92649. To insure prompt repair, be sure that you provide complete legible information and a note indicating the nature of the product issue. During the warranty period, we will repair or replace (at our sole discretion) your product at no charge. Zadro Health Solutions, Inc. shall not be liable for any incidental or consequential damages arising out of defects in this product, whether or not caused by its negligence.

Some states do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on how long an implied warranty lasts, so the above limitations or exclusions may not apply to you. This warranty gives you specific legal rights and you may also have other rights, which vary from state to state.